FCC WARNING

This device complies with part 15 of the FCC rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

NOTE: This equipment has been tested and found to comply with the limits for Class B digital devices, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determind by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and the receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

GB INTRODUCTION

Pjama Speaker is a part of an Alarm System used to treat bedwetting (enuresis nocturna) in children from the age of 5 years. The treatment can usually be carried out in the child's home. The Speaker (Audio device) connects to Pjama Sensor.

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
GB PJAMA SPEAKER

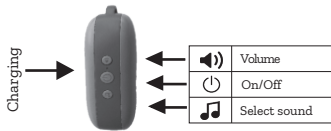INSTRUCTION

For best results: The child should be motivated and want to get dry. The parent needs to help the child wake up, especially at the beginning of treatment.
NOTE: Charge the Audio Device before use. Use the supplied USB-C cable.

Charging via:	Time- fully charged:
Computer	4 - 6 hours
AC Adapter 5V 1.0A	3 hours

During charging a red LED light flashes. When Charging is complete a red LED light will be lit constantly.

HOW TO OPERATE

1. Turn on the Speaker by holding down the  button on the Audio Device. A blue LED light starts flashing to the left of the  button and you will hear a sound that the alarm is activated. The red LED flashes when the unit needs charging.
2. Press the  button and select one of the seven different alarm sounds. The last played sound will be the sound alarming. Adjust the volume with the button .
3. Test the alarm unit. Short-circuit the Sensor by using a metallic object. The range between the Sensor and the Audio Device is approximately 10m. If your alarm does not seem to function correctly, see the instructions under 'Pairing - Pjama Sensor'.

WARRANTY


The Pjama Speaker has a one- year warranty from the buyer's purchase date of the alarm. You must be the original buyer and show evidence of your purchase. This warranty does not apply if malfunction is caused by an accident or similar incident or injury, water damage, negligence, abnormal use or failure to comply with the instructions in this manual.

CAUTION

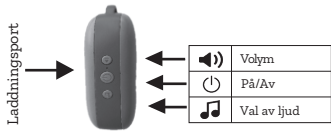- Pjama Speaker is intended for use as part of an enuresis alarm.
- Make sure that no part of the Speaker is put in the mouth for any reason. This could cause choking or other harm.
- The Pjama Speaker is a medical device that should be treated with care.
- Do not hold an alarming Audio Device within 10 cm (4 inches) of the ear. This could cause damage to the ear.
- Keep the device dry and away from direct sunlight, dust, and chemicals.

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